

# Jazzy air<sub>2</sub><sup>®</sup>

## EMPOWERING



**THE FUTURE OF MOBILITY IS HERE.**

### ***Elevate your life to the next level!***

At the touch of a switch, the Jazzy Air<sup>®</sup> 2 elevates 12" in just 11 seconds. With a sleek, more modern look, this power chair safely drives at speeds up to 4 mph while elevated, allowing users to enjoy face-to-face social engagement. With Mid-Wheel 6<sup>®</sup> Drive and Active-Trac<sup>®</sup> Suspension, the Jazzy Air<sup>®</sup> 2 offers enhanced performance and stability both indoors and out.



[pridemobility.com](http://pridemobility.com)

# Specifications

<b>Weight Capacity</b>	300 lbs.
<b>Maximum Speed<sup>1</sup></b>	Up to 4.0 mph (elevated and non-elevated position)
<b>Ground Clearance<sup>2</sup></b>	1" at motor bracket
<b>High-Back Seat to Floor Range</b>	21.25" – 33"
<b>Turning Radius<sup>2</sup></b>	20.75"
<b>Overall Size<sup>2,6</sup></b>	Length: 41.93" Width: 23.19"
<b>Drive Controller</b>	PG VR2
<b>Drive Wheels</b>	10" flat free
<b>Caster Wheels</b>	6" flat free
<b>Suspension</b>	Active-Trac <sup>®</sup> with Mid-Wheel 6 <sup>®</sup> Technology
<b>Range<sup>1,4</sup></b>	Up to 16.8 miles (U1) Up to 19.2 miles (40 Ah)
<b>Component Weights<sup>5</sup></b>	Base: 230 lbs. <sup>7</sup> w/ U1 batteries 245 lbs. <sup>7</sup> w/ 40 Ah batteries U1 Batteries: 24.5 lbs. 40 Ah Batteries: 32.5 lbs.
<b>Drivetrain</b>	Two motor, Mid-Wheel 6 <sup>®</sup> drive
<b>Brakes</b>	"Intelligent Braking" (electronic regenerative, disc park brakes)
<b>Batteries<sup>3</sup></b>	(2) 12 V, U1 (2) 40 Ah
<b>Battery Charger</b>	5 amp, off-board
<b>Warranty</b>	5 year limited (frame) 13 month limited (batteries) 13 month limited (electronics) 13 month (drive motors)

1) Varies with base model, user weight, terrain type, battery amp hour rating (Ah), battery charge, battery condition, motors, controller type, tire type, and tire condition. This specification is subject to a variance of +10%, -5%.

2) Due to manufacturing tolerances and continued product improvement, this specification can be subject to a variance of (+ or -) 3%.

3) AGM or Gel-Cell type required.

4) Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

5) Battery weight may vary based on manufacturer.

6) Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.

7) Weight includes base, seat and batteries.

NOTE: This product conforms to all applicable ANSI/RESNA, ISO 7176 Series, and EN12184 standards. All specifications subject to change without notice.



Find a compatible lift for your vehicle online:  
[PRIDEMOBILITY.COM/LIFT-FINDER](http://PRIDEMOBILITY.COM/LIFT-FINDER)

# Features

- High-back or solid seat pan, reclining sport seat. Seat features 2" of depth adjustment. Choice of 16"x16"-18", 18"x18"-20", 20"x20"-22"
- Angle, height and depth adjustable armrests
- Height adjustable headrest
- Newly re-designed lift features:
  - 2" lower seat-to-floor height enables user to sit comfortably at standard table or desk
  - Raises and lowers in 11 seconds by using the Air switch or VR2 joystick
- Additional features:
  - Drives safely at 4 mph in the elevated position
  - Integrated foot platform with rubber insert and 2" of adjustment for comfort in both the elevated and non-elevated positions
  - Bright LED marker lights offers visibility in all environments
  - Active-Trac<sup>®</sup> suspension with Mid-Wheel 6<sup>®</sup> Technology allows for maximum stability in the elevated and non-elevated positions
  - Single-sided caster forks for enhanced maneuverability
  - Under seat front storage compartment keeps your valuables accessible

# Accessories

- Alternative size lap belt: 50", 70", 80", 90"
- Cane holder
- Cell phone holder
- Cup holder
- Flag kit
- Oxygen holder
- Rear basket
- Saddle bag
- Walker holder
- Weather cover
- XLR USB Charger

# Options

- 40 amp AGM batteries
- Swing-away joystick
- Stealth<sup>®</sup> cushions

# Colors



Citrine Yellow



Sapphire Blue



Pink Topaz




Onyx Black



Ruby Red



Pearl White

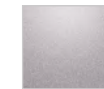

Silver

# New Matte Colors



Black Pearl



Tanzanite



Garnet Red

Authorized Pride<sup>®</sup> Provider